



KL-D.III.WMI

WICKLAND ELECTRIC OPERATING TABLE

KL-D.III.WMI

KL-D.III.WMI electric operating table is a comprehensive multi-function, which can meet the needs of various operations.



MULTIFUNCTIONAL
Multiple functions and positions



SMOOTH
Smooth operation and low volume



DURABLE
Professional quality certification

SPECIFICATIONS

- (1) Table length: $2000 \pm 100\text{mm}$
- (2) Table width: $520 \pm 20\text{mm}$
- (3) Height Adjustment: $(700-1100) \pm 50\text{mm}$
- (4) Tilt(left)/ Tilt(right): $\geq 18^\circ / \geq 18^\circ$
- (5) Trendelenburg /Rev-Trend: $\geq 25^\circ / \geq 35^\circ$
- (6) Head board upwards: $\geq 50^\circ$
- (7) Head board downwards: $\geq 90^\circ$
- (8) Back board upwards: $\geq 70^\circ$
- (9) Back board downwards: $\geq 30^\circ$
- (10) Waist board lifting: $100 \pm 20\text{mm}$
- (11) Leg board Upwards: $\geq 10^\circ$
- (12) Leg board Downwards: $\geq 90^\circ$
- (13) Leg board Outwards: $\geq 180^\circ$
- (14) Lumbar upper corner " \wedge " : $\leq 150^\circ$
- (15) Lumbar lower corner " \vee " : $\leq 110^\circ$
- (16) Horizontal and longitudinal movement: $300 \pm 30\text{mm}$
- (17) Main Voltage AC: $220\text{V} \pm 22\text{V}, 50 \pm 1\text{Hz}$
- (18) Input power VA: $400 \pm 10\%$

BASIC CONFIG

- (1) Host: 1
- (2) Screen frame: 1
- (3) Power cord: 1 pc
- (4) Arm rest: 1 pc
- (5) High-level shoulder support: 1
- (6) Shoulder support: 1 pc
- (7) Waist support: 1 pc
- (8) Leg support: 1 pc
- (9) Mattress: 1 set
- (10) Handheld Controller: 1 set
- (11) Bondage belt: 1 set



- The tabletop adopts high quality of carbon plastic plate
- The mattress is wrapped with imported artificial leather. High-density memory sponge is formed in one piece, with no seam.
- Hand held controller and panel controller are available to control various movements.
- Hydraulic system from HOERBIGER
- AC and DC can use at the same time. It has the function of auto-restoration.
- One press to " \wedge " or " \vee ".
- Imported sealing ring
- The base and rails on both sides of the table are 304 stainless steel.
- Electric braking device, convenient for shifting the table.

MAIN FEATURE

KL-D.III.WMI



PRODUCT DESCRIPTION



Safe and convenient operation

- Handheld remote control and panel controller in two ways.
- Self-locking function to improve clinical safety
- With built-in casters, electric brake device, convenient operation table shift

High-quality composite material

- The base, table support frame and side rails are made of durable and easy-to-clean 304 stainless steel
- The bed surface and backboard adopt imported carbon-plastic panels, which can meet the C-arm body inspection and X-ray filming.



Imported humanized mattress

- The mattress is wrapped with imported artificial leather, and the high-density memory foam is molded at one time without gaps
- Soft cushion ensures that the patient is evenly stressed
- Anti-static, waterproof, detachable, easy to clean



MULTI-SCENARIO APPLICATION



NEUROLOGY




OBSTETRICS



INTERNAL MEDICINE



ENT

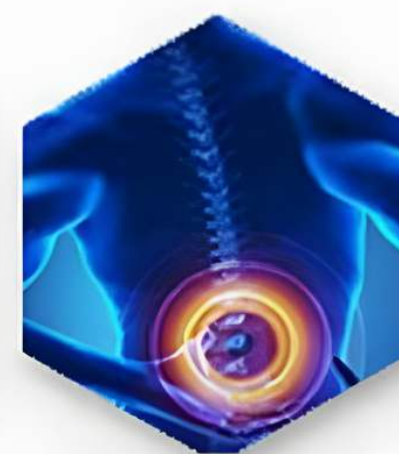

ORTHOPEDICS



PT. Wickland Medika Imantara

Medical, Aesthetics & Laboratory Supplier

Jl. Palem Sememi Raya No 12 – 12A Kel. Sememi,
Kec. Benowo Surabaya 60198 | 031-33100412