



PT. Wickland Medika Imantara

Medical, Aesthetics & Laboratory Supplier

Jl. Palem Sememi Raya No 12 - 12A Kel. Sememi,
Kec. Benowo Surabaya 60198 | 031-33100412



BR-E11 Series
FLEXIBLE VIDEO
BRONCHOSCOPE

▼ Specifications

BR-E11PF

Insertion Tube	2.8mm
Distal End	3.1mm
Working Channel	1.2mm
Working Length	610mm
Field of View	120°
Depth of Field	3-100mm
Bending(Up/Down)	210°/130°
Rotation(Left/Right)	120°/120°
Application	Pediatric

BR-E11F

Insertion Tube	3.7mm
Distal End	3.9mm
Working Channel	1.7mm
Working Length	610mm
Field of View	120°
Depth of Field	3-100mm
Bending(Up/Down)	210°/130°
Rotation(Left/Right)	120°/120°
Application	Pediatric/Pneumology

BR-E11PT

Insertion Tube	3.9mm
Distal End	4.1mm
Working Channel	2.0mm
Working Length	610mm
Field of View	120°
Depth of Field	3-100mm
Bending(Up/Down)	210°/130°
Rotation(Left/Right)	120°/120°
Application	Pediatric

BR-E11H

Insertion Tube	4.8mm
Distal End	5.2mm
Working Channel	2.0mm
Working Length	610mm
Field of View	120°
Depth of Field	3-100mm
Bending(Up/Down)	210°/130°
Rotation(Left/Right)	120°/120°
Application	Pneumology

BR-E11HT

Insertion Tube	5.8mm
Distal End	6.0mm
Working Channel	2.8mm
Working Length	610mm
Field of View	120°
Depth of Field	3-100mm
Bending(Up/Down)	180°/130°
Rotation(Left/Right)	120°/120°
Application	Pneumology

BR-E11TB

Insertion Tube	5.8mm
Distal End	6.2mm
Working Channel	3.0mm
Working Length	610mm
Field of View	120°
Depth of Field	3-100mm
Bending(Up/Down)	180°/130°
Rotation(Left/Right)	120°/120°
Application	Pneumology



Suction Valve Retainer

Unique design prevents suction valve from falling off.



Multi-Function Buttons

- ① Zoom in / out
- ② Capture / record
- ③ Freeze / release



Rotation

Improve operation flexibility and accuracy.



Passive Bending

Angle Lock

Bending angle locking function for easy diagnosis and treatment.



Wide Bending Angle

smaller bending radius.

▼HV-3101+BR-E11 Series

▼HV-3030+BR-E11 Series

Image Processor HV-3101

- ▶ 10.1" Touch screen
- ▶ 1280*800 HD resolution
- ▶ HDMI+AV output
- ▶ Rechargeable battery supporting 4 working hours
- ▶ Portable design



Image Processor HV-3030

High-pixel image quality, multi-function integration, and highly compatible image processing platform.

1. Automatic Gain Control:




Adjusts image brightness for optimal visual quality.

2. Contour Enhancement:



Boosts image clarity to suit doctors' varying needs for detail.

3. Light Metering Mode:

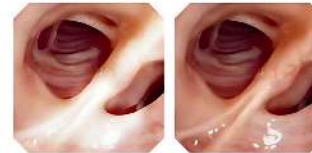

Adjusts exposure to eliminate glare, enhancing image quality.

4. Customizable Buttons:



Enables quick and efficient operation for enhanced diagnosis and treatment.



- ▶ 4.0" independent touch screen
- ▶ Video standard: 1920 × 1080
- ▶ Output interface: DVI * 2 SDI * 2 CVBS * 2
- ▶ Memory: 64GB SD card

